

Operating Instruction


- 1.Insert plug in power outlet.
- 2.Turn on the Power Switch:
- 3.Choose functions on Control Panel.
 - Massage**: Push the button up for kneading massage.
 - Infrared**: Push the button up for infrared heat.
 - Kneading and heat function can work together or separately.
- 4.When treatment is completed, turn off the power button on Power Switch.
- 5.Unplug from outlet.

Specifications:

Name: Massage Cushion

Input: ☐ AC220-240V~, 50/60Hz
☐ AC100-120V~, 50/60Hz

Power: 55W

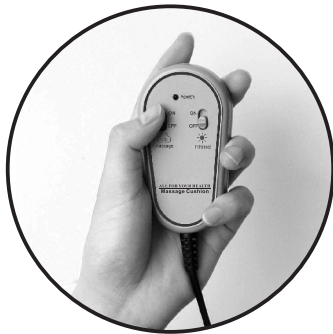

Control Panel

Massage Cushion

It is 6 kneading massage balls inside the cushion with 2 infrared rays devices.



Instruction Manual



WARNING: For proper use, the user must read and understand the operator's manual before using this product.

SAVE THIS MANUAL FOR FUTURE REFERENCE

Important Safety Instructions

When using an electrical appliance, basic precautions should always be followed including the following:

DANGER

To reduce the risk of electric shock:

Always unplug this appliance from the electrical outlet right now after using or before cleaning.

WARNING

To reduce the risk of burns, fire, electric shock, or injury :

- 1) An appliance should never be left unattended when plugged in.
- 2) Do not operate under blanket or pillow. Excessively heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this massager is used by, on or near children, invalids or disabled persons.
- 4) This massager is only for whom needed as described in this manual. Do not use attachments that not recommended by the manufacturer.
- 5) Never operating the massager if it has a damaged cord/plug, not working properly or dropped into water. Return the massager to a service centre for test and repair.
- 6) Do not carry this massager by supply cord or use it as a handle.
- 7) If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- 8) Never operating the massager with the air openings blocked. Keep the air openings free of lint, hair, etc.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operating where aerosol (spray) products are being used or where oxygen is being administered.
- 12) Keep the cord away from heated surfaces.
- 13) To disconnect, turn all controllers to the off position, then remove plug from outlet.

SAVE THESE INSTRUCTIONS

Consistent and acute massage treatment has been proven to be an effective way to improve general well-being. Massage helps to

- ◆ Reduce muscular tension and associated discomforts.
- ◆ Enhance tissue elasticity and flexibility.
- ◆ Relieve muscle aches and stiffness.
- ◆ Reduce muscle soreness after exercise.
- ◆ Improve athletic performance.
- ◆ Reduce stress and anxiety.

- ◆ Improve circulation and oxygenation.
- ◆ Strengthen immune system.
- ◆ Release toxins.



Please stop using this product and consult your doctor as soon as possible if you have one of the following symptoms:

- Severe pain.
- Swollen or painful joint.
- Wound, bruise or broken bone.

Do not use this product as a substitute for medical treatment!

Product Functions

- The 6 finger-like massage nodes perform kneading massage technique.
- To massage your neck, up or lower back along both sides of the spine.
- Massage nodes rotate and bring vital nutrients to spinal area.
- 2 lenses emit infrared rays that accelerate the kneading action which in turn allows heat to penetrate deeply into muscle layers in order to improve circulation and ease tension.
- With this massage cushion, you can initiate your massage experience virtually anywhere: sofa, bed, recliner, chair, or even on the floor.

A short period of massage treatment will lift up your spirit!

